

MDX Series



Model
MDX-101
Brix / Refractive Index
Refractometer
USER GUIDE



Included Parts



1. MDX-Series Digital Refractometer
2. Removable Rubber Grip
3. Removable Prism Cover
4. Disposable Pipet
5. AAA, 1.5V Battery (1ea.)
6. Hardshell Storage Case
7. Screwdriver (not shown)

Instrument Use

Power On and Calibration

One press of the power key will power on the instrument (figure 6). In order to maximize battery life the refractometer is designed to shut itself off after 1 minute of inactivity.

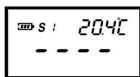


figure 6

A zero calibration should be performed daily to achieve optimal accuracy.

- Clean and dry the sample well. Place a 4-5 drops of clean, distilled water in the sample well.
- Close the prism cover to prevent stray light from affecting the calibration.
- Press and hold the ZERO key for two seconds until "CAL" is displayed, then release (figure 7).
- "CAL" will flash repeatedly for 10 seconds. Press ZERO once again while "CAL" flashes to perform the calibration. During this time it is crucial to not open the prism cover.
- If ZERO is not pressed again then the calibration will not be performed.
- The display will indicate "END" when finished calibrating (figure 8).



figure 7



figure 8



figure 9



figure 10

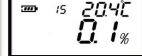


figure 11

- At this point the display will change back to the default screen showing the active scale and instrument temperature (figure 9). If the instrument does not read distilled water with the values shown below then clean the prism and recalibrate with fresh distilled water:

Distilled Water
0.0% Brix ±0.2%
1.3330RI ±0.0003RI

Scale Selection

Repeatedly press and release the SCALE key to toggle through all available scales. The display will indicate which scale position is currently active- S01, S02, S03,... (figure 9).

You can find descriptions and ranges for the available scales on the bottom of the instrument, the storage case label, and this user guide.

Taking Measurements

Single Readings:

- Ensure that the sample well is clean and dry. Place 4-5 drops of your liquid sample into the well and close the prism cover.
- Press the READ key once and the sample reading will display on the screen (figure 10). Additional readings can be made by pressing the READ key again if desired.
- Clean and dry the sample well when finished. The prism must be cleaned and dried before testing another sample.

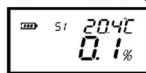


figure 10

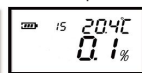


figure 11

Introduction

The VEE GEE MDX-Series portable digital refractometers are a modern design utilizing advanced opto-electronic components to provide years of fast, accurate liquid testing with high reproducibility.

The MDX-Series models are designed to be simple to use with little training necessary, regardless though, it's highly recommended that users read this User Guide in its entirety before using the instrument for the first time.

At any point if you have any questions please contact our support team at:

Warranty

VEE GEE digital refractometers are warranted to be free from defects in material and workmanship for a period of one year from the date of purchase. During this period the VEE GEE Service Center will, at their option and without charge, either repair or replace any part found to be defective in materials and workmanship.

All warranty work shall be performed by the VEE GEE Service Center. Contact your dealer or the VEE GEE Service Center for troubleshooting the issue and to receive a Return Authorization, if necessary, for the return of your instrument for repair. The party returning the product must prepay all postage, shipping, transportation, packaging, duties & taxes, and delivery costs to the VEE GEE Service Center.

This warranty is subject to the following limitations and will not apply if:


- 1) There is lack of proof of purchase date and place of purchase. The warranty is not assignable or transferable.
- 2) The damage is due to normal wear (including the prism), misuse, abuse, negligence, or any other cause not due to manufacturing of the product.
- 3) The serial number is altered or obliterated; or unauthorized repair or replacement of parts by any party other than the VEE GEE Service Center.

This warranty expressly excludes transportation damage and readjustment. In no case shall VEE GEE Scientific LLC be liable to the Buyer or any person for any special, indirect, incidental, or consequential damage whether claims are based in contract or otherwise with respect to or arising out of product furnished hereunder. For goods manufactured by any third party, VEE GEE Scientific's liability under warranty is limited to the terms of the warranty by the supplier of the goods.

3


VEE GEE MDX-Series

Instrument Components

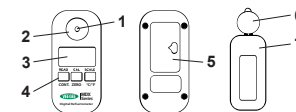

1. Prism Surface
2. Sample Well
3. LCD Display
4. Control Keys
5. Battery Compartment Cover
6. Removable Prism Cover
7. Removable Rubber Grip

Instrument Display

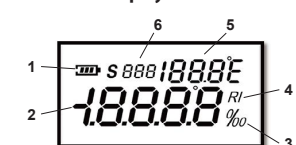


1. Battery Level Indicator
2. Reading Display
3. Percent or Per mille Scale Indicator
4. Refractive Index Scale Indicator
5. Temperature Display (°C or °F)
6. Active Scale Indicator / Continuous Read Timer

4

| Battery Level | Display Indicator |
|---------------|-------------------|
| 100-80% | |
| 80-50% | |
| 50-20% | |
| <20% | |

Expected battery life: 1000+ reads


Battery Compartment

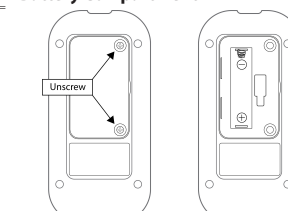

- Unscrew the two screws to open the battery compartment.
- Place AAA, 1.5V (1ea.) battery into the battery slot. Pay careful attention to ensure proper battery polarity when inserting.

5

Control Keys

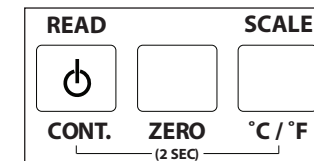

READ | Power On | CONT.

- Press once and release to power on.
- Press once and release while powered on to take reading.
- Press and hold for 2 seconds while powered on to initiate continuous read loop function.

ZERO

- Press and hold for 2 seconds then press once again within 10 seconds to perform zero calibration procedure.

SCALE | °C / °F

- Repeatedly press and release to toggle through all available scales.
- Press and hold for 2 seconds to toggle between °C and °F temperature display.

6

Error Codes



figure 12

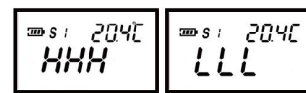


figure 13



figure 14

| Error Code | Description |
|-------------------------|--|
| A01 | Calibration Temperature Error: Operating temperature must be 0-40°C. (figure 12) |
| A02 | Calibration Sample Error: Use clean distilled water only for zero calibration. |
| A03 | Hardware Failure: Contact our Technical Support department. |
| HHH / LLL (sample) | Sample High/Low Value Error: Sample above or below measurable range. (figure 13) |
| HHH / LLL (temperature) | Temperature High/Low Value Error: Environment above or below acceptable range. (figure 14) |

10

| Model | MDX-101 |
|--------------------|-------------------|
| Catalog Number | 48101 |
| Scale 1 Range | 0.0 - 50.0% Brix |
| Scale 1 Resolution | 0.1% Brix |
| Scale 1 Accuracy | ±0.2% Brix |
| Scale 2 Range | 1.3330 - 1.4200nD |
| Scale 2 Resolution | 0.0001nD |
| Scale 2 Accuracy | ±0.0003nD |

11

| Model | MDX-101 |
|----------------------|------------------------------------|
| ATC Range | 0-40°C 32-104°F |
| IP Rating | IP67: detecting unit IP65: body |
| Wavelength | 589.3nm LED |
| Sample Cell | Stainless Steel |
| Calibration Solution | Distilled Water |
| Min. Sample Volume | 0.3mL |
| Battery | 1 X AAA, 1.5V |
| Dimensions | 121 x 58 x 25mm |
| Weight | 90g |
| Warranty | 1-Yr Parts & Labor |



8

9