


REVISION RECORD			
LTR	DESCRIPTION	DATE	BY
A	RELEASED FOR PRODUCTION TYPO: ADDED 3" DIMENSION	5,1,17 5,11,17	28218 bs/ bs/

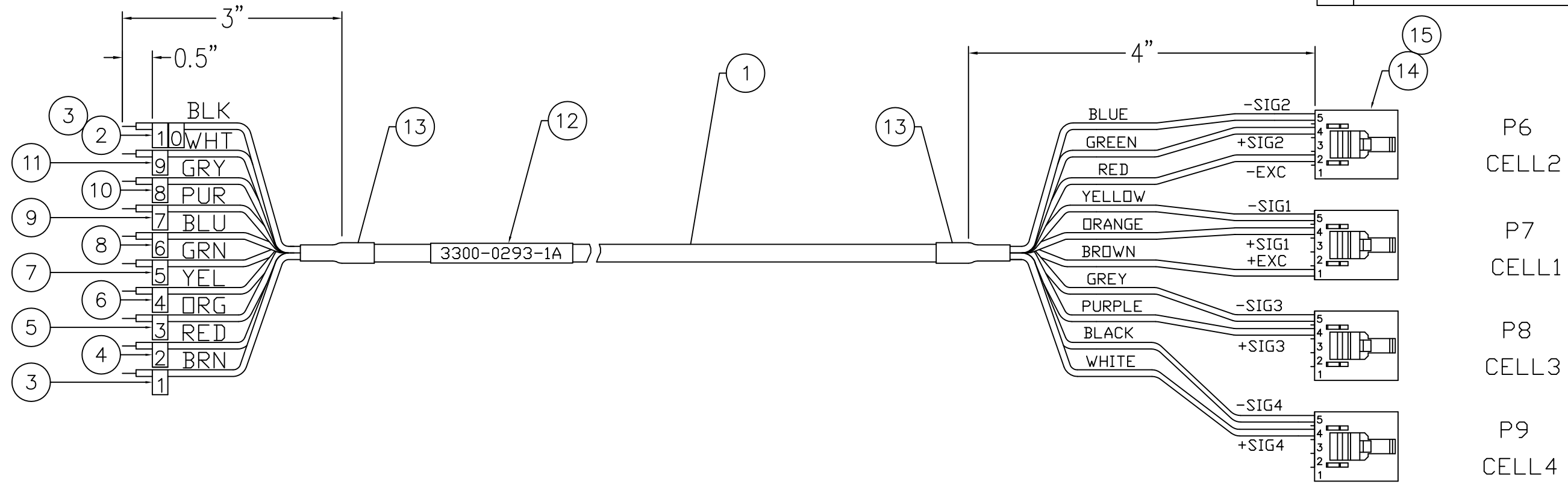

NOTE:
1. STRIP EACH WIRE .25 AND TIN.

10	BLACK	-SIG4	P9-5
9	WHITE	+SIG4	P9-4
8	GREY	-SIG3	P8-5
7	PURPLE	+SIG3	P8-4
6	BLUE	-SIG2	P6-5
5	GREEN	+SIG2	P6-4
4	YELLOW	-SIG1	P7-5
3	ORANGE	+SIG1	P7-4
2	RED	-EXC	P6-2
1	BROWN	+EXC	P7-1

10	15	6610-5182	TERMINAL, FEMALE	
4	14	6610-5181	CONNECTOR, 5 PIN	FEMALE
.083'	13	6980-0114	HEAT SHRINK	
1	12	6650-0078	LABEL	
.028	11	6980-0149	WIRE MARKER 9	
.028	10	6980-0148	WIRE MARKER 8	
.028	9	6980-0147	WIRE MARKER 7	
.028	8	6980-0146	WIRE MARKER 6	
.028	7	6980-0145	WIRE MARKER 5	
.028	6	6980-0144	WIRE MARKER 4	
.028	5	6980-0143	WIRE MARKER 3	
.028	4	6980-0142	WIRE MARKER 2	
.056	3	6980-0141	WIRE MARKER 1	
.028	2	6980-0140	WIRE MARKER 0	
10.0'	1	6980-0099	CABLE	

QTY.	QTY.	QTY.	ITEM	PART/DWG. NUMBER	DESCRIPTION	MATERIAL OR SOURCE	WT. EACH
			WT.	PARTS/MATERIAL LIST			

UNLESS OTHERWISE SPECIFIED TOLERANCE ON DIMENSIONS ARE:
 ANGLES ± 1/2°
 INTEGERS/FRACTIONS ± 1/16 IN.
 DECIMALS (.X) ± .03 IN.
 DECIMALS (.XX) ± .01 IN.
 DECIMALS (.XXX) ± .005 IN.
 NOTE: WHOLE NUMBERS MUST BE WRITTEN XX.000 TO INVOKE DECIMAL TOLERANCES.

Cardinal
 Cardinal Scale Manufacturing Co.
 WEBB CITY, MISSOURI, USA 64870

CARDINAL DETECTO FULLER WEIGHING SYSTEMS

TITLE: DIALYSIS LOAD CELL CABLE

SCALE: N.T.S. DO NOT SCALE DRAWING NEXT ASSEMBLY

DR. bs DATE 5/1/17

CH. DATE

MODEL ID- SHEET 1 OF 1 DWG. NO. 3300-0293-1A REV. A

THIS DRAWING IS THE PROPERTY OF THE CARDINAL SCALE MFG. CO. AND SHALL NOT BE USED NOR REPRODUCED EITHER WHOLLY OR IN PART EXCEPT WITH WRITTEN AUTHORIZATION. ALL RIGHTS OF DESIGN OR INVENTION ARE RESERVED.