

BLICKMAN



Case Carts

Compact, stable, secure transport


Includes:

- Closed Case Carts
- Open Case Carts
- Shelving Options
- Accessories

THE BLICKMAN DIFFERENCE

- Operating room surgical staff can transport fully loaded case carts throughout the facility with the doors closed for increased protection and cart cleaning is a simple and easy task.
- We offer options that meet a variety of different surgery needs including neurosurgery, general surgery, outpatient surgery, labor and delivery, and more. We also offer space saving solutions that take up less floor space while being able to handle more extensive procedures.
- Our broad offering of case carts remains the most durable and highest functioning on the market. We are the only manufacturer in our class to offer a lifetime limited warranty. Our advanced engineering capabilities enable us to customize case cart designs to satisfy virtually any requirement your medical facility may have.
- The doors of the closed case carts swing a full 270° and a positive spring acting clip at the cart ends capture doors to stay in the open position during cart washing and while in use in the operating room.

Key Product Features

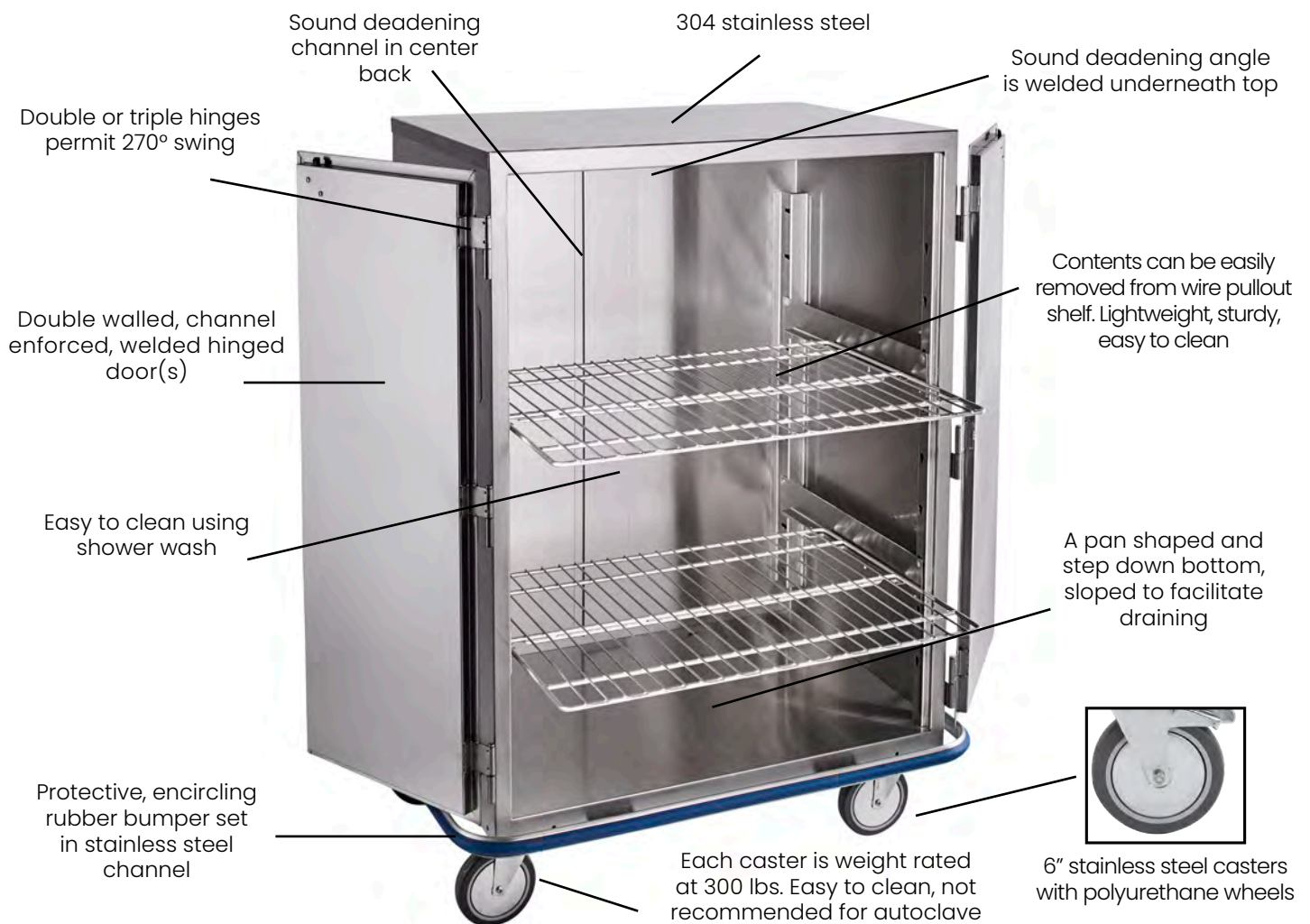

Case Carts

Maxi to Mini

Description: Our white maxi, multipurpose, and mini case carts differ in width, they are all 40 1/2" tall and 29" deep.. They each come standard with one wire shelf. The maxi case cart comes standard with a heavy duty roll-out shelf, and the multipurpose and mini case carts come with a lightweight pull-out shelf.

Maxi Cart

52"W

Ideal for procedures that require a larger volume of supplies such as neurosurgery and orthopedic surgery

Internal Dimensions:
42 7/8"W x 29 3/8"H x 25"D

Multipurpose Cart

42" W

Ideal size for most routine surgical procedures, whether or not your facility uses a container system

Internal Dimensions:
32 7/8"W x 29 5/8"H x 25"D

Mini Cart

29 5/8" W

Designed for labor and delivery and outpatient surgery, where a smaller cart is preferred

Internal Dimensions:
20 1/2"W x 29 5/8"H x 25"D

Push Handle



Name

Maxi

Multipurpose

Mini

Solid Doors

Model #	CCC1-19	CCC2-19	CCC3-19
Item #	2293331000	2293332000	2293333000
Floor Space	10 1/2 ft ²	8 1/2 ft ²	6 ft ²
Shipping Weight	266 lbs.	225 lbs.	168 lbs.

Glass Doors

Model #	CCC1-19G	CCC2-19G	CCC3-19G
Item #	22933G1000	22933G2000	22933G3000
Floor Space	10.47 ft ²	8.46 ft ²	5.97 ft ²
Shipping Weight	260 lbs.	219 lbs.	165 lbs.



Case Carts

Space Saving

Description: Our space saver case carts are designed to give a greater capacity while using the same floor space as the mini and multipurpose case carts. Their height and additional shelving provide additional storage capacity. Instead of 40 1/2" tall, they are 55 1/2" tall, and instead of one shelf, they come standard with two lightweight pull-out wire shelves.

Space Saving Cart

29 5/8" W

This cart can handle most routine surgical procedures.

Ultra Space Saving Cart

42" W

This cart can handle more extensive procedures.

Recessed Handle



Internal Dimensions:

20 1/2"W x 44 5/8"H x 25"D



Internal Dimensions:

32 7/8"W x 44 5/8"H x 25"D



Name

Space Saving

Ultra Space Saving

Solid Doors		
Model #	CCC4-19	CCC5-19
Item #	2293334000	2293335000
Floor Space	6 ft ²	8 1/2 ft ²
Shipping Weight	208 lbs.	276 lbs.
Glass Doors		
Model #	CCC4-19G	CCC5-19G
Item #	229333G4000	229333G5000
Floor Space	6 ft ²	8 1/2 ft ²
Shipping Weight	198 lbs.	262 lbs.

Case Carts

Over The Road

Description: As healthcare processes change, Blickman works with customers to design the products they need. With growing demand for off-site sterile processing, Blickman reengineered its case cart lineup by adding a new series of carts dedicated to meet its customers' new and more demanding transport requirements.



**Ultra Space Saver - OTR5 is shown here with optional vertical side handles.*



Name	Item #	Model #	Dimensions	Floor Space	Shipping Weight
Maxi	2327331030	OTR1	42H x 53W x 29D	10 1/2 ft ²	266 lbs.
Multipurpose	2327332030	OTR2	42H x 43W x 29D	8 1/2 ft ²	225 lbs.
Mini	2327333030	OTR3	42H x 30 5/8W x 29D	6 ft ²	168 lbs.
Space Saver	2327334030	OTR4	56 3/4H x 30 5/8W x 29D	6 ft ²	208 lbs.
Ultra Space Saver	2327335030	OTR5	56 3/4H x 43W x 29D	8 1/2 ft ²	276 lbs.

Case Carts

Open

Description: Our open case cart designs are easy to clean, dry, load, and maneuver. The full encircling bumper and solid bottom prevents cross contamination. This design is ideal for carts that go directly to a sterile core, bypassing non-sterile transport.

Overall Dimensions:

Width: 42" x Height: 42" x Depth: 29"



Name	Solid Shelves	Wire Shelves
Model #	OCC3	OCC4
Item #	2410396000	2410496000
Shipping Weight	219 lbs.	141 lbs.

Case Carts

Accessories

ACCESSORIES

Clean/Dirty Indicators



US Design Patent #: US D883,378 S
US Utility Patent#: US 10,643,505 B1



US Design Patent #: US D883,379 S
US Utility Patent#: US 10,643,505 B1

Part #	Description	Shipping Weight
5278628KIT	Angle Model	2 lbs.
5200045KIT	Dial Model	2 lbs.
5278618KIT	Pod Model	2 lbs.



Card Holder

We offer a 8.5" x 11" card holder. Custom sizes and placement options available. Stainless steel frame withstands prolonged use and cutouts allow for drainage.



Tote Boxes

Sturdy, long lasting tote boxes are ideal for organizing and storing small supplies that may slip through the spaces of wire shelves.

Part #	Description	Shipping Weight
5206701185	Card Holder	N/A - Factory Installed
9854194720	Tote Boxes	5 lbs.



Case Carts

Construction Specifications

Construction Specifications

Universal Dimensions	Closed Case Carts: 29" deep, varying widths. Standard Height: 40 1/2"; Space saving height: 55 1/2" Open Case Carts: 42" x 42" x 29"
Materials	304 stainless steel, non-magnetic
Gauge	Closed Case Carts: Body & Doors: 18GA; Top & Chassis: 16GA Open Case Carts: Top & Bottom: 16GA; Shelf: 18GA Bumper Channel: 14GA
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	Four 6" stainless steel swivel casters with polyurethane wheels, weight rated at 300lbs (one with directional lock)
Notes	Found on next page

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs.

Case Carts

Construction Specifications

Notes Continued

CLOSED CASE CARTS:

- **Body:** One piece top and one piece body. Sound deadening angle welded under the top. Sound deadening channel on the center back. A pan shaped, step down bottom with a slope that facilitates draining. Horizontal intersections filled with heat resistant silicone. The body shell is bolted to the chassis. A tubular push handle is provided at one end of the cart.

- **Doors:** Double walled, channel enforced, (double or triple) hinged doors with paddle latch pulls on leading edges. Doors less than 36" high swing on two hinges, while doors greater than 36" high swing on three hinges. All doors permit full 270° swing. In fully opened position, positive spring acting clips at the cart ends capture doors for cart washing or operating room use.

- **Shelves:** Shelves slide on adjustable stainless steel wires, suspended from keyhole punched standards welded to the cart interior.

OPEN CASE CARTS:

- **Frame:** U-shaped tubular framing formed around 3" radii or tubular corner posts stabilized and capped by a pan shaped top welded to posts at contact points. Posts are fitted with casters and either rolling bumpers or encircling rubber bumpers.

- **Tops, bottoms, fixed shelves:** Fabricated with turned down edges, either hemmed or flanged, with corners notched to sit over or snugly against the inside contour of the upright members. Notched corners are fully welded to the uprights the full height of the facing edges as well as around the notch.

- **Adjustable shelves:** Solid rods welded between the tubular uprights at each end, stabilizing the framing members while providing varying levels of shelf adjustment.
