

# HEAT SEALED POWDER COAT STEEL 10 SCOPE CAPACITY CABINET, KEY LOCK



**SureDry**  
SCOPE CABINETS with

**Dri-Scope Aid®**  
**TRACKER**

Scopes not included

**SC8030TD-T2310D**

*Also available  
in a Stainless Steel Finish*

## FEATURES

- Positive pressure fully redundant fan system
- Unit features two intake fans and one exhaust fan
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Holds five scopes with secure adjustable mounting system
- Replaceable HEPA Filter system without removal of scopes
- Medical grade electrical system
- Tempered glass double doors
- Multiple locking options available including badge entry
- Easy to clean removable drip tray
- Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada

## CABINET DIMENSIONS

- 93" H x 30" W x 24" D

## COLORS

Available in more than twenty color options.



**Includes Dri-Scope Aid®, industry leading internal lumen drying system**

Please refer to Dri-Scope Aid® specification sheet for more information

**HARLOFF**  
INNOVATING STORAGE FOR CARE

Engineered &  
Manufactured  
in the U.S.



The ultimate in Endoscope drying systems is now available! SureDry Scope Cabinets feature an external HEPA filtered positive pressure drying system. Dri-Scope Aid® is the industry leading internal lumen drying system. This combination provides facilities with the assurance that they will be able to meet all current and anticipated future guidelines for external and internal drying of endoscopes, providing their patients the best possible assurance of infection prevention.

# DRI-SCOPE INTERNAL LUMEN DRYING SYSTEM W/TRACKER



**SureDry**  
SCOPE CABINETS with

**Dri-Scope Aid®**  
**TRACKER**

## TRACKER FEATURES & BENEFITS


- Digital display maintains a list of scopes by drying status, with displays of scope processing times
- Intelligent automated scheduler continuously processes scopes as long as they remain checked into the system
- Includes default or customizable settings and drying options unique to each endoscope
- Scope tracker maintains records of scope drying events, user access, check-ins and check-outs
- Color touch screen access by authorized personnel is completed by passcode, with adjustable session timeout settings and options for passcode retrieval if forgotten
- Various adjustable timing options controlled through the color touchscreen including drying, pre-drying and venting time settings, hang time limits, endoscope removal/load delay timing, etc
- At a glance, display shows scope inventory and reprocessing status including the time until scope expiration requiring further reprocessing
- Slot display order settings displayed on the status screen per user preferences
- Configurable replacement times for tubing and pump filter, with overdue warnings/indicator lights
- Downloadable export files available for endoscope management, user lists and more

The Dri-Scope Aid® Cabinet supplies HEPA-filtered air to all internal channels to aid in drying endoscopes with an automated drying cycle for improved turn-around time, and a continuous venting cycle for overnight storage. The Dri-Scope Aid® system in the SureDry Scope Cabinets provides a complete drying system for your facility.



### Customizable Hang Time Monitor

Adjust the number of days your scopes hang. You will be alerted when a scope is about to expire

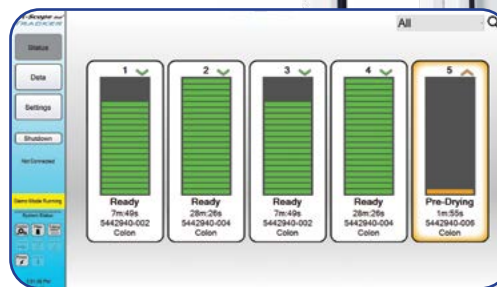

**HARLOFF**  
INNOVATING STORAGE FOR CARE

Engineered &  
Manufactured  
in the U.S.



**Optional dimmable LED lighting also available!**



# 10 Reasons to Choose **HARLOFF** for Your Scope Cabinets



**Industry Best Positive Pressure HEPA Filtered Drying System.** Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

1

*A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines*

## Proudly Made in the USA



Harloff products are proudly **Designed, Engineered and Manufactured in the United States** at our Colorado Springs, CO facility.

6




## LED Lighting Option

**Brighten your cabinet** and more easily inspect scopes with fully dimmable optional overhead LED lighting which triggers on/off when you open and close the door.

2

## Easy Access to HEPA Filters

Change HEPA filters **without opening the cabinet.** Save time and ensure no cycle interruptions.

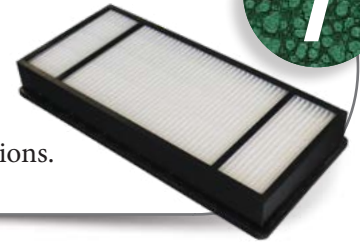

7

## Disposable Connectors

**Save time and ensure compliance!** Use single-use connectors for all internal channels of the endoscope. Part of the Dri-Scope Aid® package.



3

## Keyed or Keyless Locks

**Options include** key locks, basic electronic push button locks, advanced keypads and proximity readers.



8

## Hooks and Holders for Specialized Scopes

Major brands to include **Fuji, Olympus, Pentax, Karl Storz and others.**



4

## Finish Options and Colors

Stainless Steel, Heat Sealed Powder Coat (High Gloss or Textured Low Gloss inside & out)

9

## 12 Year Warranty

**Our Manufacturer's Warranty is Unmatched** in the industry. We build our cabinets to last and we stand behind them.

5

## Mobile Scope Cabinet Kit

Convert any size of scope cabinet into a mobile version. Includes 4 3" casters and steel mounting plate.



10