



model 8440 BABY SCALES



Detecto
Model 8440
Baby Scale

Providing accurate and safe weighing of infants and toddlers, Detecto's model 8440 digital scale is ideal for home health and medical clinics. The 8440's weighing tray may be removed for standing toddler weighing and for easier transport. The 1-in / 25-mm high LCD offers pounds and kilograms readouts plus features a built-in clock. The scale's Weight Difference key provides the weight difference between the current weight on the scale and a saved weight in memory along with its time and date.

Detecto
Model
8440-CASE
Carrying
Case






Large, easy-to-read display with weight recall



The tray collapses for easier transport and to allow for convenient standing toddler weighing.

- Accuracy to 1/2 ounce
- Storage of previous weight memory recall
- Calculation of change in weight
- Low battery indicator
- Optional carrying case available (8440-CASE)

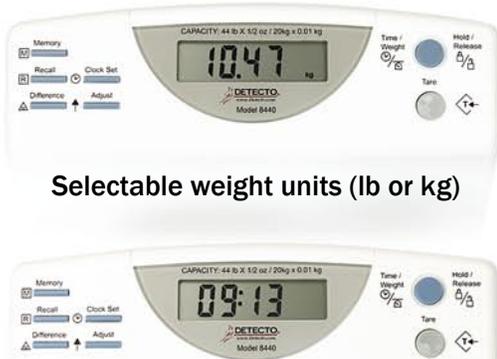

Selectable weight units (lb or kg)

Date and time displays automatically in sleep mode or by pressing the Time/Weight button during weighing

8440 BABY SCALE SPECIFICATIONS

CAPACITY:	44 lb x 1/2 oz / 20 kg x 10 g
DISPLAY:	1 in / 25 mm high LCD
TRAY DIMENSIONS:	20.6 in W x 10.5 in D x 4 in H / 523.9 mm W x 266.7 mm D x 101.6 mm H
SCALE DIMENSIONS:	9.8 in W x 11.3 in D x 7.5 in H / 247.7 mm W x 285.8 mm D x 190.5 mm H
SHIPPING WEIGHT:	8 lb / 3.6 kg
POWER:	(1) 9V Alkaline battery (not included) and (1) 3V CR2032 lithium back-up battery (included) for memory
CLOCK FUNCTION:	Time and date
MEMORY FUNCTION:	Memory for comparison of current and stored weight

Detecto reserves the right to improve, enhance, or modify features and specifications without prior notice.



Sold By:



www.universalmedicalinc.com