

THE 553KL


DIGITAL PEDIATRIC SCALE WITH EXTRA-WIDE TRAY

The 553KL pediatric scale was specifically created for nurseries, delivery rooms, and pediatric practices, making it an ideal choice for newborns, infants, and toddlers. With a maximum capacity of 44 pounds and 20 kilograms, this scale offers providers and patients a range of benefits, including safety, comfort, and flexibility. The scale ensures a secure platform, keeps patients calm for accurate readings, and accommodates infants lying down to toddlers sitting up.

The 553KL pediatric scale offers a built-in measuring tape, allowing users to measure weight and height simultaneously. It also provides the flexibility of choosing between an optional digital or mechanical height rod. For easy portability and storage, there is an optional carrying case.

With its user-friendly digital interface, the 553KL scale includes features such as Zero, Tare, Hold/Release, and Auto-off, making it convenient and easy to use. Its motion detection technology ensures accurate weight measurements even when the patient moves.

The 553KL includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

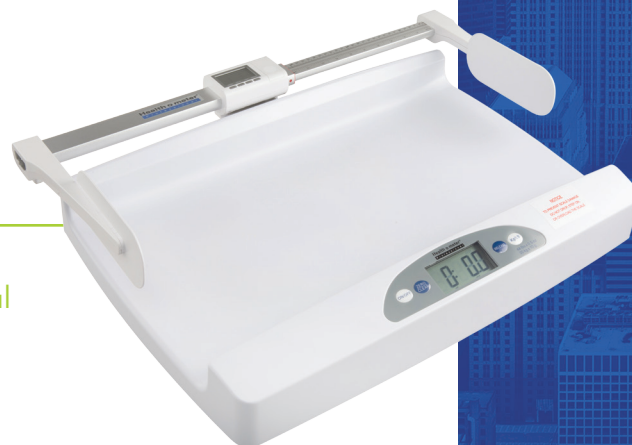

553KL-EHR

Scale with attachable digital height rod

THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 553KL is your choice for safety, comfort, and reliability.

120+ LOCAL
FIELD SPECIALISTS

100%
CONNECTIVITY

#1 TRUSTED
BRAND

100 YEARS
OF SERVICE

Health o meter®
Professional

Health o meter® Professional 553KL

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

44 lb / 20 kg

Resolution

0.5 oz/ 10g

EMR Connectivity

Via USB, optional Pelstar® wireless technology

Display

1" / 25 mm, LCD display

Power Source

6 AA batteries included, 100-240V adapter (ADPT30) not included

Functions

LB / KG conversion

LB / KG lock out

Zero

Tare

Hold / release

Auto zero

Auto off

Motion detection

Tray Size (w x d x h)

24 7/8" x 13" x 1 3/8"
(632 mm x 330 mm x 35mm)

Product Footprint

24 7/8" x 16 1/2" x 3 7/8"
(632 mm x 419 mm x 98 mm)

Product Weight

10 lb / 5 kg

Measuring Tape

0" - 24" / 0 cm - 58 cm

Portability

553CASE - optional carrying case

2210CART- optional rolling cart

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

MODEL	DESCRIPTION
553KL	Digital Pediatric Scale with Extra-Wide Tray, LB & KG
553KG	Digital Pediatric Scale with Extra-Wide Tray, KG only
553KL-HR	553KL with attachable mechanical height rod
553KL-EHR	553KL with attachable digital height rod
553HR	Mechanical height rod for 553 series of scales
553EHR	Digital height rod for 553 series of scales
SS-553KL	ScaleSurance® 2 year extended warranty for 553 series*
2210CART	Rolling scale/ medical device cart for 553 series of scales
553CASE	Carrying case for 553 series of scales
ADPT30	Optional adapter for 553 series of scales
PLUGSET30	Plug adapters for use outside North America
PELSTARONE-P	Wireless thermal printer



553CASE

Carrying Case for 553 Series of Scales



2210CART

Medical Device Cart for 553 Series of Scales

Health o meter®
Professional



*Not available in all countries.