



P-2040

blanket warming cabinet specifications

- 15.4 ft³ (0.44m³) interior capacity
7.7 ft³ (0.2m³) per chamber
- Adjustable temperature setting up to
98°F - 200°F (37°C - 93°C)
- Two compartments with individual
temperature control
- Digital display
- Epoxy-coated blanket support assembly
- Stainless steel exterior

Specifications:


Dual-cabinet design with 20 gauge stainless steel exterior casing and doors. Window doors allow observation of inventory with the doors closed. The cabinet is warmed using Radiant Heat Technology low-heat-density electrothermal cable array, mounted against the walls, eliminating the need for a heat circulating fan. Each compartment is equipped with one white, epoxy-coated blanket support assembly and one center shelf. The cabinet is furnished with one set of 5" (127mm) heavy duty casters, two with locking brake.

Control:

Each compartment is controlled by one electronic control which consists of a 4 digit L.E.D. display, on/off key, heat indicator L.E.D., increase and decrease keys, built-in lockout feature and a series of prompt sequence indicators. Each control has an adjustable temperature range of 98° to 200°F (37° to 93°C). Controls can easily be set to operate in Celsius or Fahrenheit. Timer key allows the user to program the control to automatically turn on and turn off once during a 24 hour period at a selected time. In the event of a power failure, the cabinet will remember its programming and begin to operate as before. A warming shut-off system, separate from the electronic control, prevents overheating.

Special Information:

- Right-hand door swing is standard. Left-hand door is available as a special order.
- Pass-through design is available.

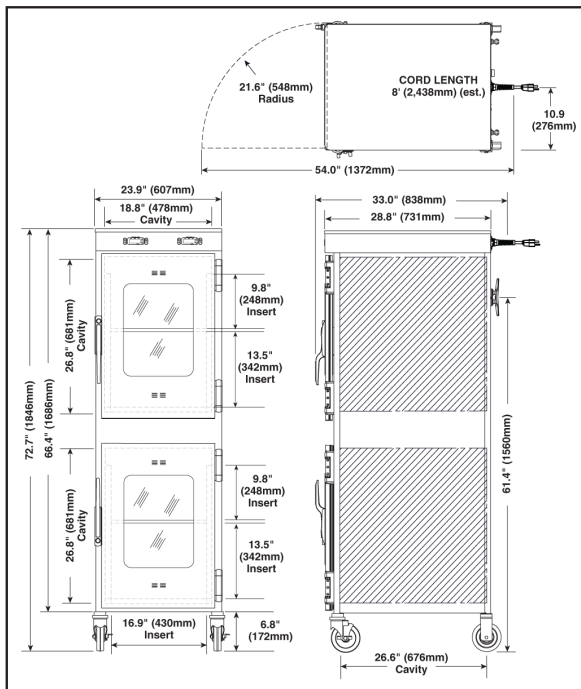

Capacity & Weights	
Interior Capacity	7.7 cubic feet (0.2m ³)
Net Weight	412 lb. (187 kg) est.
Ship Weight	502 lb. (228 kg) est.
Electrical	
120 V.A.C. 50/60 Hz, 1 ph 2000 Watts, 16.0 Amps	NEMA 5 - 20 P 20A - 125V Plug
230 V.A.C. 50/60 Hz, 1 ph 1800 Watts, 7.8 Amps	BS 1363 Plug CEE 7/7 (UK only) 220-230V
Options & Accessories	
Cylinder Lock for Door Handle	LK-22567
Full-Perimeter Rubber Bumper Assembly	5001763
Leg Kit, 6" (152mm), Flanged	44093
Window Door Assembly, Left Hand (each)	5001255