



Task & Examination Light

Gleamer®

Burton's all-purpose halogen light.

If your facility is searching for an all-purpose, rugged, compact halogen light, Burton's Gleamer® is your solution. The Gleamer® was designed to meet healthcare facilities' most demanding specifications for high-intensity, cool operation, compact design and value. Offering 45,000 lux, a sturdy gooseneck, plus optional 16" (41 cm) extension arm, and superior quality components, the Gleamer® provides high-intensity, maximum flexibility and extraordinary reliability. The Gleamer® also features a unique venting system, energy-absorbing infrared exchanger and reflector coatings, which filter out infrared for a nearly heat-free beam in a cool, yet compact housing.

Aesthetics, advanced optics and ergonomics ensure uncompromising performance and comfort.



The Right Light


Features

- Easy bulb replacement
- Cool housing temperature and cool light beam; light can be in continuous use for hours without causing discomfort to physician or patient
- Convenient positioning with ergonomic positioning ridges
- High-quality, space-saving 4-caster base – safe and flexible portable movement with locking feature in 2 of 4 casters



Summary of technical data	
Illuminance:	45,000 lux (4181 fc) at 18" (46 cm)
Color temperature:	3200 K (natural coloring)
CRI (Color Rendering Index):	91
Type of bulb:	35-watt halogen bulb; 2000-hour bulb life
Power requirements:	120V model standard; 230 / 240V model available
Certifications/Approvals:	IEC/UL 60601-1 / IEC 60601-2-41 / IEC/EN 60601-1-2; CAN/CSA C22.2 601-1 M90 certified CE marked
Product warranty:	5-year limited warranty
Origin:	Assembled in the USA

Note: The optical data are nominal values based on measurements calculated according to UL/IEC performance standards at 120V

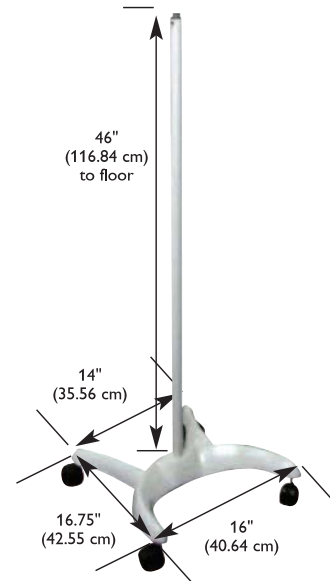

Sample Model No: GLW32FL = Gleamer® Standard Spot, 230 / 240V, Floorstand

Gleamer®	Voltage	Model
GL = Standard Spot GLW = Wide Beam Spot	30 = 120V 32 = 230 / 240V	FL = Floorstand W = Extension Arm and Universal Wall Mount

Accessories

- #0009600PK Gleamer® Replacement Bulbs, 4/Box
- #0009608PK Gleamer® Wide Beam Spot Replacement Bulbs, 4/Box